Portable massage table

Bärbar massagebänk / Bærbart massasjebord /
Tragbare Massagetisch / Kannettava hierontapöytä

ENG

SE

NO

DE

FI

Thank you for choosing to purchase a product from Rusta!

Read through the entire manual before installation and use!

Portable massage table

Read the user manual carefully to make sure that the product is assembled, used and maintained correctly, as described in the instructions in this user manual. Keep the user manual for future reference.

SAFETY INSTRUCTIONS

- Check the table periodically. Make sure all screws are secure and the cable has not deteriorated or frayed. Make sure there are no cracks in the wood structure.
- Do not overload the table by exceeding the recommended maximum working weight (205kg) - the combined weight of the client and maximum amount of downward pressure being applied by the therapist. Doing so may result in serious injury and/or damage.
- Do not use alcohol-based cleaners, as it could cause damage to the vinyl materials. Instead, use warm soapy water and wipe dry.
- Do not expose the table to extreme (hot or cold) temperature or moisture.
- To ensure stability and safety, always position your client in the center of the table to distribute the weight evenly. Also make sure that, when in use, all four legs are on a solid ground and at the same level.


WARNING!


If these instructions are not followed, an accident may lead to serious personal injury or property damage.

SET-UP INSTRUCTIONS

1. Place the table on the side, release the locks, partially open the table, and remove the accessories.

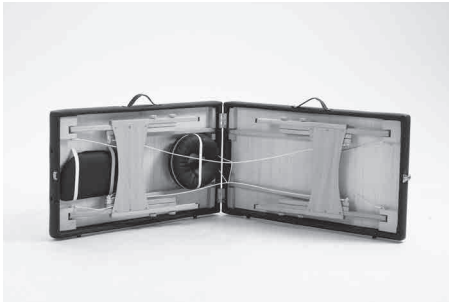
2. Carefully extend the legs, making sure the cables are not tangled.

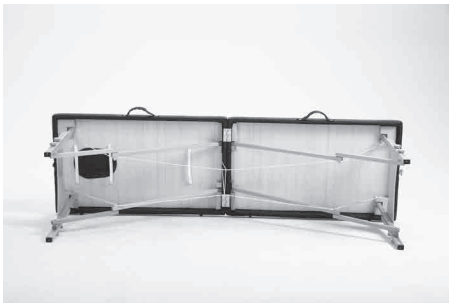3. Turn the table upright by grasping the edge of the table in the middle and pulling the table up. Adjust the table to make sure the surface flattens completely.

4. Adjust the height of the table by removing both knobs on each table leg. Then place the knobs back securely in the desired position.

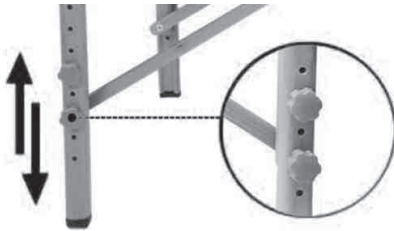HEADREST AND ARM SHELF SET-UP INSTRUCTIONS


1. Insert the headrest tubular into the headrest holes in the table and lock the position using the locking lever of the headrest.

2. Place the remaining single arm shelf strap through the middle part of headrest.
3. Thread the remaining single strap through the buckle.

4. Adjust the arm shelf to the desired horizontal position and tighten all the straps accordingly.
5. Attach the headrest pillow to the headrest with the provided Velcro fastener.

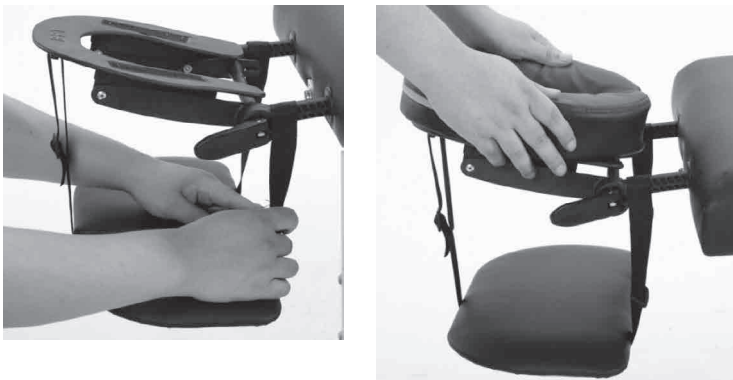ADJUSTING THE HEADREST

1. Position the locking lever so that it extends out horizontally with parts on top.

2. Adjusting the headrest as desired and then move the levers downwards to lock into place.

CARE AND CLEANING

1. When not in use, keep the table in the carrying case and store it in a cool dry place, avoiding direct sunlight.
2. For daily cleaning, remove daily dirt from the table with mild soap and warm water and dry it with a soft, line-free cloth.
3. For stain removal, use a mild, non-abrasive cleaner only. Spray a cloth with cleaner and wipe the stain. Rinse thoroughly with warm water, and dry it immediately.
4. The hardwood and is finished with lacquer, thus it requires no special care. Keep the table dry and avoid exposure to wet or humid environments at all time.

TECHNICAL DATA

Dimensions	L 185 x W 55 x H 61–88 cm
-------------------	---------------------------

RIGHT OF COMPLAINT

By law the product must be returned to the place of purchase along with the original receipt in the event of a complaint. The user is liable for any damage to the product that is caused by not using the product for its intended purpose or by not following this user manual correctly. The right of complaint will not apply in these instances.

INFORMATION ABOUT HAZARDOUS WASTE

Used products should be handled separately from household waste. They should be taken for recycling in accordance with local rules for waste management.

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