# Neck & shoulder massager

Massageapparat nacke och axlar / Nakke- og skuldermassasje /  
Nacken- & Schultermassagegerät / Hierontalaite niskalle ja hartioille

ENG

SE

NO

DE

FI

Thank you for choosing to purchase a product from Rusta!

Read through the entire manual before installation and use!

## Neck & shoulder massager

Read the user manual carefully to make sure that the product is assembled, used and maintained correctly, as described in the instructions in this user manual. Keep the user manual for future reference.

## SAFETY INSTRUCTIONS

- Do not allow any foreign matter to enter the clearance between the massage head and the casing. Do not allow the product to come into contact with water so as to avoid any possible accidents that may result from a short circuit.
- If the mains cable or plug is damaged, do not disassemble and attempt to repair it by yourself. Please contact your distributor.
- Ensure that the mains power is suitable for this product.
- Please disconnect the product from the mains supply after use or before cleaning to avoid damage to the equipment or personal injury.
- Use only the AC adaptor provided.
- Do not use any spare parts and attachments other than those recommended.
- Do not use this product if its casing or mains cable is damaged.
- Do not use this product if you are under the influence of alcohol or feeling unwell.
- Do not use this product with too much force so as to avoid injury.
- Do not keep this product on for more than 15 minutes so as to avoid overheating.
- Do not use this product within one hour of eating.
- Please use this product carefully according to your doctor's advice if you have a neck or back injury or if you have had an operation recently.
- Do not plug or unplug this product when your hands are wet.
- Do not pull the cable when unplugging. Do not plug or unplug carelessly.
- This appliance is not intended for use by persons (including children)

with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be done by children unless they are older than 8 and supervised. Keep the appliance and its cord out of reach of children less than 8 years.
- Do not damage the cable or change the circuit of the product.
- Please cease using this product in case of a power cut in order to avoid any body injury upon sudden power resumption.
- Please stop using this product and consult your local supplier if you notice anything abnormal with the product during use.
- If the supply cord is damaged, it should be replaced by supplier, retailer, it's service agent or a similarly qualified person in order to avoid a hazard.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- Please do not pull controller connecting line, in cause that power cord fracture or plug with bad contact. When pulling out power plug, please hold the plug pulled back, rather than directly draw the power cord.
- When the whole machine does not working, please immediately stop using it to pull out power plug. This product is prohibited in the following: Bathroom and kitchen high temperature lights, damp, high voltage and the source of fire and water.
- Do not force it with fierce pressure or put body pressure on the massage mat. Don't touch with hard objects; don't weigh or tread it. Never interact with the zipper.

## **WORKING ENVIRONMENT**

- Do not use this product in a highly humid environment such as a bathroom.
- Do not use this product immediately after there has been a sharp change in the environmental temperature.
- Do not use this product in an environment where there are lots of dust or corrosive gases.
- Do not use this product when it is covered by a blanket or pillow, etc.

## **PEOPLE THAT SHOULD NOT USE THIS PRODUCT**


- People with contagious skin conditions.
- Patients with osteoporosis.
- Patients with heart diseases and with electronic medical instruments like a heart pacemaker.
- People in a high fever.
- Women that are pregnant or menstruating.
- Those who are injured or suffering from any body surface disease.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- Those who to need rest according to a doctor's advice or those who do not feel well.
- Do not use this product when you are wet.
- Do not use this product on a pet.

### **WARNING!**

**If these instructions are not followed, an accident may lead to serious personal injury or property damage.**

## PRODUCT FEATURES

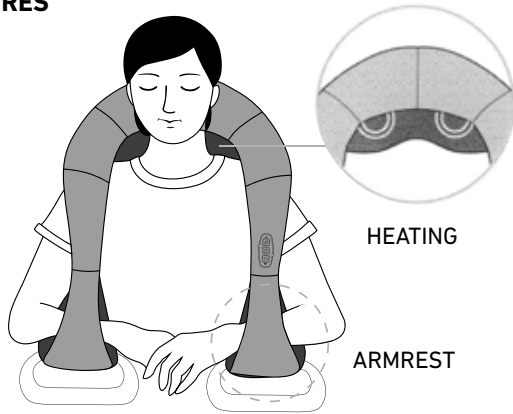- The new concept of kneading massage, matches the perfect fit of the neck and shoulder.
- The massage head is specially designed according to the neck and shoulder to give you the most comfortable experience.
- U shaped, arm in arm design, strong or gentle type, strength is decided by the user.
- Clockwise and anti-clockwise kneading automatic conversion system.
- Kneading direction manually chosen function.
- Warm function controll separately, promote the blood circulation effectively.
- Polyurethane high-grade material holster.
- Equipped with a overheating safety protection device.

## METHOD AND FUNCTION OF BUTTONS

### HT (HEATER)

You can choose to stop heating function when it is working, press one time, stop heating function, press again, heating function work again.

### ON/OFF: ON or OFF and Kneading Direction Control

Press to Start working, automatically changing kneading direction every minutes, press again to set only one direction, then another direction by the third time. Hold for a few seconds to stop the massage.

### SPEED CONTROL

3 levels of speed, the default is high speed, press one time is low speed, press again middle speed, then high speed.


### ROTATION

The product will automatically rotate clockwise or counterclockwise every minute. Push the Rotation button to stop the auto-rotate mode. Then you can change the direction by manually push the Rotation button.


Waist massage

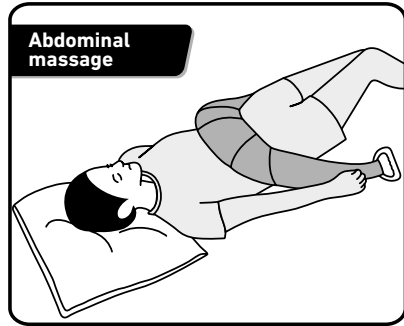Abdominal massage



Leg massage

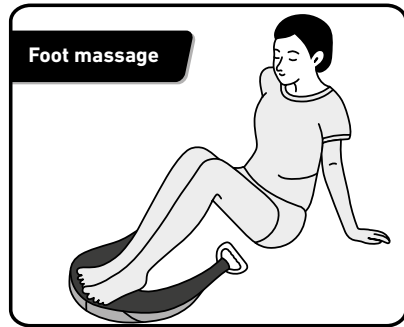Foot massage

### LOCATION OF PARTS AND CONTROLS

1. Heat button
2. ON/OFF and Kneading direction button
3. Speed control button
4. Rotation button
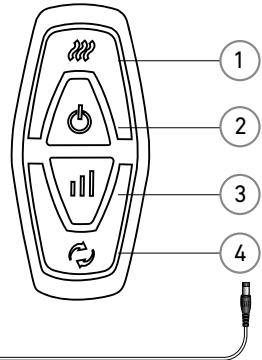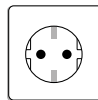
### POWER

1. Plug adaptor in to the massage appliance.
2. The power plug insert into the socket.

### CLEANING AND MAINTENANCE

1. After using the products, please unplug it from the mains before you start cleaning.
2. Cleaning is performed with a cloth for wiping of dirt from the product.
3. Please do not use corrosive cleaner, such as benzene or thinner to clean this product. Use only water with a cloth.
4. Do not use a wet cloth, to clean the control panel, use only a dry cloth.
5. Products when not in use for long time, please wipe clean, place in well ventilated cases and dry place, in order to avoid the product be affected with damp and affect use.

## TROUBLESHOOTING

- This chapter sums up the most frequently encountered problems in using this product.
- If you have problems not mentioned in the list, you should try the following steps.
- Turning off the power of the product, unplug the power plug from the socket.
- Let it cool down for one 1 hour.
- Insert the plug into the power socket and open the switch.
- If the problem persists, please contact our customer service center.

<b>Trouble shooting</b>	<b>The reason</b>	<b>Management</b>
The product can't work.	The socket is unplugged.	Plug the socket safely.
	Power button is not activated.	Activate the power button.
	This product has just finished massage 20 minutes and automatically switch machine to standby mode.	ON/OFF. If you want to massage again, please click ON/OFF again.
	Products placed on the overload (for example, someone standing on the products, or have a weight placed on top of the product).	Unplug the product power plug, remove the excess weight. Please insert product power and will be open again.
Product overheated. The power supply cord damaged.	You might have exceeded the amount of massage cycles in a row.	Remove the power plug, let it cool down for 1 hour, if still does not work, contact the seller.
The thermostat heating disconnect protection.	Internal short circuit.	Return to the seller.

**TECHNICAL DATA**

<b>Adapter Input</b>	100–240V, 50/60Hz
<b>Adapter Output</b>	DC12V, 2A, 24W

**WARRANTY**

Please keep the proof of purchase (receipt) before claiming warranty.

1. This product is manufactured under strict quality control and inspection processes.
2. When this product goes out of order under the normal use of customer, we will repair it free of charge during the warranty period.
3. Please show us the warranty card when you ask for repair service.
4. Please keep the warranty card securely as it is not reissued.
5. Shape or circuit of the product are subject to change without any notice to improve the quality of product.
6. The damages (defaults) prescribed below are NOT to be covered by warranty.

Users faults by lack of care.

Unauthorised electric currency connection defaults.

Faults by the users own intention of dismantle or repairs.

Damages caused by natural disaster (fire, flood, earthquake, lighting, etc.).

Replacement to new parts.

## **RIGHT OF COMPLAINT**

By law the product must be returned to the place of purchase along with the original receipt in the event of a complaint. The user is liable for any damage to the product that is caused by not using the product for its intended purpose or by not following this user manual correctly. The right of complaint will not apply in these instances.

### **INFORMATION ABOUT HAZARDOUS WASTE**

Used products should be handled separately from household waste. They should be taken for recycling in accordance with local rules for waste management.

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