Lumbar massager

Massagekudde / Massasjepute /
Massagegerät / Hierontatyyny

ENG

SE

NO

DE

FI

Thank you for choosing to purchase a product from Rusta!

Read through the entire manual before installation and use!



Lumbar massager

Read the user manual carefully to make sure that the product is assembled, used and maintained correctly, as described in the instructions in this user manual. Keep the user manual for future reference.

SAFETY PRECAUTIONS

- Do not use your massager in an excessively humid or dusty environment, as it may cause electrical shock or machine trouble. Do not use the massager near hot surfaces or areas.
- Keep the massager away from children.
- The following users must consult their doctor before use:
 - Those who are under medical treatment.
 - Those who suffers from hypertension, heart disease, brain disease, tumours and/or thrombosis.
 - Pregnant women.
 - Those who are suffering from osteoporosis.
 - Those who have pacemakers or other implanted medical devices.
- Do not operate the massager without carefully reading this manual.
- Do not allow children to play around the massager. It may cause injury.
- Please leave the massager to cool down for 20 minutes after each use.
- Stop using if you are feeling unwell. Please stop using if water spills over the massagers, if any parts are damaged, if it starts to smell or smoke during operating
- Do not use while sleeping.
- Do not connect the massager with any other power supply or electrical equipment at the same time.
- Keep away from water and any other liquids.
- Do not place the massage balls directly on your hands, elbows, knees or chest.
- The massager may get hot on the surface. Those sensitive to heat shall pay special attention to it

DESCRIPTION

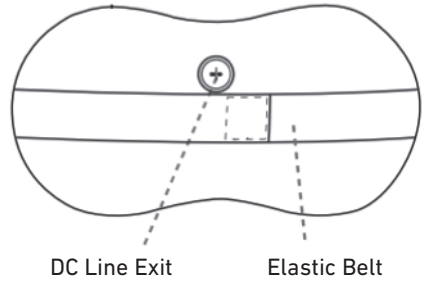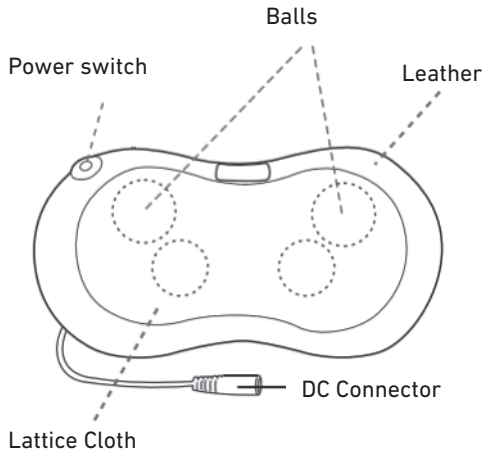DC 12V Power Adaptor

Car Adaptor Charges

WARNING!

If these instructions are not followed, an accident may lead to serious personal injury or property damage.

INSTRUCTIONS

- Insert the power plug into proper wall outlet.
- Press the ON/OFF button to start Massage Pillow. Press again to switch on heating function.
- The massage balls operate clockwise and anticlockwise a minute at a time before rotating.
- Press ON/OFF button again to switch off the heating function. Press it again if you want to start the heating function.
- Hold the ON/OFF button in 3 seconds to turn the massager off completely.

SPECIFICATIONS

Output voltage	DC 12V
Power	24W
Net weight	1.35 kg
Gross weight	1.75 kg

RIGHT OF COMPLAINT

By law the product must be returned to the place of purchase along with the original receipt in the event of a complaint. The user is liable for any damage to the product that is caused by not using the product for its intended purpose or by not following this user manual correctly. The right of complaint will not apply in these instances.

INFORMATION ABOUT HAZARDOUS WASTE

Used products should be handled separately from household waste. They should be taken for recycling in accordance with local rules for waste management.

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